

ACQUITY UPLC HSS COLUMNS

CONTENTS

I. GETTING STARTED

- a. Column Connectors
- b. Column Installation
- c. Column Equilibration
- d. eCord Installation
- e. Initial Column Efficiency Determination
- f. VanGuard™ Pre-Columns

II. COLUMN USE

- a. Sample Preparation
- b. pH Range
- c. Solvents
- d. Pressure
- e. Temperature

III. COLUMN CLEANING, REGENERATING AND STORAGE

- a. Cleaning and Regeneration
- b. Storage

IV. INTRODUCING ECORD INTELLIGENT CHIP TECHNOLOGY

- a. Introduction
- b. Installation
- c. Manufacturing Information
- d. Column Use Information

V. ADDITIONAL INFORMATION

- a. Tips for Maximizing ACQUITY UPLC HSS Column Lifetimes
- b. Recommended Flow Rates and Backpressures for Reversed-Phase ACQUITY UPLC HSS Columns

Thank you for choosing a Waters ACQUITY UPLC® HSS column. The ACQUITY UPLC HSS packing materials were designed specifically for use with the Waters ACQUITY UPLC system and are manufactured in a cGMP, ISO 9001:2000 certified plant using ultra pure reagents. Each batch of ACQUITY UPLC HSS material is tested chromatographically with acidic, basic and neutral analytes and the results are held to narrow specification ranges to assure excellent, reproducible performance. Every column is individually tested and a Performance Chromatogram and Certificate of Batch Analysis are provided on the eCord™ intelligent chip.

ACQUITY UPLC columns were designed and tested specifically for use on the ACQUITY UPLC system. ACQUITY UPLC columns will exhibit maximum chromatographic performance and benefits ONLY when used on the holistically-designed ACQUITY UPLC system since the ACQUITY UPLC system and column were created and designed to operate together. For these reasons, Waters cannot support the use of ACQUITY UPLC columns on any system other than an ACQUITY UPLC system.



I. GETTING STARTED

Each ACQUITY UPLC HSS column comes with a Certificate of Analysis and a Performance Test Chromatogram embedded within the eCord intelligent chip. The Certificate of Analysis is specific to each batch of packing material contained in the ACQUITY UPLC HSS column and includes the gel batch number, analysis of unbonded particles, analysis of bonded particles, and chromatographic results and conditions. The Performance Test Chromatogram is specific to each individual column and contains such information as: gel batch number, column serial number, USP plate count, USP tailing factor, capacity factor, and chromatographic conditions. These data should be stored for future reference.

a. Column Connectors

The ACQUITY UPLC system utilizes tubing and gold plated compression screws which have been designed to meet stringent tolerance levels and to minimize extra column volumes.

Optimized column inlet tubing (Part Number 430001084) is supplied with the ACQUITY UPLC system. The inject valve end of the tubing is clearly marked with a blue shrink tube marker. Insert the opposite end of the tubing into the ACQUITY UPLC column and tighten the compression fitting using two 5/16-inch wrenches.

For information on the correct column outlet tubing, please refer to the relevant detector section in the ACQUITY UPLC System Operator's Guide (Part Number 71500082502).

b. Column Installation

Note: The flow rates given in the procedure below are for a typical 2.1 mm i.d. by 50 mm length 1.8 μm column. Scale the flow rate up or down accordingly based upon the flow rate and pressure guide provided in Section V (Additional Information).

1. Purge the pumping system of any buffer-containing mobile phases and connect the inlet end of the column to the injector outlet.
2. Flush column with 100% organic mobile phase (methanol or acetonitrile) by setting the pump flow rate to 0.1 mL/min and increase the flow rate to 0.5 mL/min over 5 minutes.
3. When the mobile phase is flowing freely from the column outlet, stop the flow and attach the column outlet to the detector. This prevents entry of air into the detection system

and gives more rapid baseline equilibration.

4. Gradually increase the flow rate as described in step 2.
5. Once a steady backpressure and baseline have been achieved, proceed to the next section.

Note: If mobile phase additives are present in low concentrations (e.g., ion-pairing reagents), 100 to 200 column volumes may be required for complete equilibration. In addition, mobile phases that contain formate (e.g., ammonium formate, formic acid, etc.) may also require longer initial column equilibration times.

c. Column Equilibration

ACQUITY UPLC HSS columns are shipped in 100% acetonitrile. It is important to ensure mobile phase compatibility before changing to a different mobile phase system. Equilibrate the column with a minimum of 10 column volumes of the mobile phase to be used (refer to Table 1 for a list of column volumes). The column may be considered thermally equilibrated once a constant backpressure is achieved.

Column Length (mm)	Internal Diameter		
	1.0 mm	2.1 mm	3.0 mm
30	—	0.1	0.2
50	0.04	0.2	0.4
100	0.08	0.4	0.8
150	0.12	0.5	1.0

To avoid precipitating mobile phase buffers on your column or in your system, flush the column with five column volumes of a water/organic solvent mixture, using the same or lower solvent content as in the desired buffered mobile phase. (For example, flush the column and system with 60% methanol in water prior to introducing 60% methanol/40% buffer mobile phase).

d. eCord Installation

The eCord button should be attached to the side of the column heater module. The eCord button is magnetized and does not require specific orientation.

e. Initial Column Efficiency Determination

1. Perform an efficiency test on the column before using it. This test may consist of:
 - a) An analyte test mixture that is commonly used in your laboratory, and/or
 - b) An analyte mixture as found on the “Performance Test Chromatogram” which accompanied your column.


Note: If b) is performed, the isocratic efficiencies measured in your laboratory may be less than those given on the Waters “Performance Test Chromatogram”. This is normal. The Waters isocratic column testing systems have been modified in order to achieve extremely low system volumes. This presents a more challenging test of how well the column was packed. This guarantees the highest quality packed column. These special testing systems have been modified to such an extent that they are not commercially viable and have limited method flexibility other than isocratic column testing.

2. Determine the number of theoretical plates (N) and use this value for periodic comparisons.
3. Repeat the test at predetermined intervals to track column performance over time. Slight variations may be obtained on two different UPLC® systems due to the quality of the connections, operating environment, system electronics, reagent quality, column condition and operator technique.

f. VanGuard Pre-Columns

VanGuard Pre-columns are 2.1 mm ID x 5 mm length guard column devices designed specifically for use in the ACQUITY UPLC system. VanGuard Pre-columns are packed with the same UPLC chemistries and frits as our 2.1 mm ID UPLC columns. VanGuard Pre-columns are designed to be attached directly to the inlet side of an ACQUITY UPLC column.

Note: In order to ensure void-free and leak-free connections, the VanGuard Pre-column is shipped with the collet and ferrule NOT permanently attached. Care must be taken when removing the O-ring that holds these two pieces on the pre-column tubing.

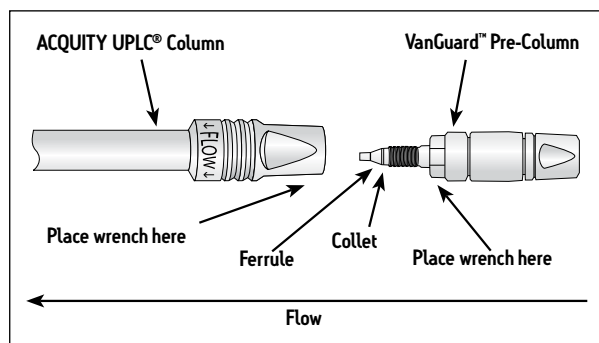

INSTALLATION INSTRUCTIONS

1. Remove the VanGuard Pre-column from its box and shipping tube and remove plastic plug.
2. Orient the pre-column so that male end is facing up and carefully remove rubber O-ring that holds collet and ferrule in place during shipping (collet and ferrule are not yet permanently attached).
3. Orient the ACQUITY UPLC column perpendicular to the work surface so that column inlet is on the bottom (column outlet on top).
4. From below, insert the VanGuard Pre-column into the ACQUITY UPLC column inlet and hand-tighten (collet and ferrule are not yet permanently attached).
5. While pushing the VanGuard Pre-Column into the column inlet, turn assembled column and pre-column 180° so that pre-column is now on top.
6. Tighten with two 5/16” wrenches placed onto the ACQUITY UPLC column flats and the VanGuard Pre-column hex nut (male end) as shown above.
7. Tighten 1/4 turn to set collet and ferrule.
8. Check that the ferrule is set by loosening the connection and inspecting the ferrule depth. A properly set ferrule depth will resemble other connections in the ACQUITY UPLC system.
9. Reattach pre-column, apply mobile phase flow and inspect for leaks.

Table 2: Buffer Recommendations for Using ACQUITY UPLC HSS Columns

Additive/Buffer	pKa	Buffer Range	Volatility (± 1 pH unit)	Used for Mass Spec	Comments
TFA	0.3		Volatile	Yes	Ion pair additive, can suppress MS signal, used in the 0.02-0.1% range.
Acetic Acid	4.76		Volatile	Yes	Maximum buffering obtained when used with ammonium acetate salt. Used in 0.1-1.0% range.
Formic Acid	3.75		Volatile	Yes	Maximum buffering obtained when used with ammonium formate salt. Used in 0.1-1.0% range.
Acetate (NH ₄ CH ₂ COOH)	4.76	3.76 – 5.76	Volatile	Yes	Used in the 1-10 mM range. Note that sodium or potassium salts are not volatile.
Formate (NH ₄ COOH)	3.75	2.75 – 4.75	Volatile	Yes	Used in the 1-10 mM range. Note that sodium or potassium salts are not volatile.
Phosphate 1	2.15	1.15 – 3.15	Non-volatile	No	Traditional low pH buffer, good UV transparency.
Phosphate 2	7.2	6.20 – 8.20	Non-volatile	No	Above pH 7, reduce temperature/concentration and use a guard column to maximize lifetime.

II. COLUMN USE

To ensure the continued high performance of ACQUITY UPLC HSS columns, follow these guidelines:

a. Sample Preparation

- Sample impurities often contribute to column contamination. One option to avoid this is to use Waters Oasis® solid-phase extraction cartridges/columns or Sep-Pak® cartridges of the appropriate chemistry to cleanup the sample before analysis. For more information, visit www.waters.com/sampleprep
- It is preferable to prepare the sample in the operating mobile phase or a mobile phase that is weaker than the mobile phase for the best peak shape and sensitivity.
- If the sample is not dissolved in the mobile phase, ensure that the sample, solvent and mobile phases are miscible in order to avoid sample and/or buffer precipitation.

Filter sample with 0.2 μ m membranes to remove particulates. If the sample is dissolved in a solvent that contains an organic modifier (e.g., acetonitrile, methanol, etc.) ensure that the membrane material does not dissolve in the solvent. Contact the membrane manufacturer with solvent compatibility questions. Alternatively, centrifugation for 20 minutes at 8000 rpm, followed by the transfer of the supernatant liquid to an appropriate vial, could be considered.

b. pH Range

The recommended operating pH range for ACQUITY UPLC HSS columns is 2 to 8 for the T3, C₁₈ SB, Cyano and PFP chemistries,

and 1 to 8 for the HSS C₁₈ chemistry. A listing of commonly used buffers and additives is given in Table 2. Additionally, the column lifetime will vary depending upon the operating temperature, the type and concentration of buffer used.

Note: Working at the extremes of pH, temperature and/or pressure will result in shorter column lifetimes.

c. Solvents

To maintain maximum column performance, use high quality chromatography grade solvents. Filter all aqueous buffers prior to use through a 0.2 μ m filter. Solvents containing suspended particulate materials will generally clog the outside surface of the inlet distribution frit of the column. This will result in higher operating pressure and poorer performance. See Section V for more information.

d. Pressure

ACQUITY UPLC HSS Columns can tolerate operating pressures up to 18000 psi (1241 bar or 124 MPa).

Note: Working at the extremes of pressure, pH and/or temperature will result in shorter column lifetimes.

e. Temperature

Temperatures between 20–45 °C are recommended for operating ACQUITY UPLC HSS columns in order to enhance selectivity, lower solvent viscosity and increase mass transfer rates. When operating at high pH, lower operating temperatures are recommended for longer column lifetime. Working at high temperatures (e.g. > 40 °C) may also result in shorter column lifetimes.

Note: Working at the extremes of temperature, pressure and/or pH will result in shorter column lifetimes.

III. COLUMN CLEANING, REGENERATING AND STORAGE

a. Cleaning and Regeneration

Changes in peak shape, peak splitting, shoulders on the peak, shifts in retention, change in resolution or increasing backpressure may indicate contamination of the column. Flushing with a neat organic solvent, taking care not to precipitate buffers, is usually sufficient to remove the contaminant. If the flushing procedure does not solve the problem, purge the column using the following cleaning and regeneration procedures.

Use the cleaning routine that matches the properties of the samples and/or what you believe is contaminating the column (see Table 3). Flush columns with 20 column volumes of solvent. Increasing column temperature increases cleaning efficiency. If the column performance is poor after regenerating and cleaning, call your local Waters office for additional support.

Table 3. Reversed-Phase Column Cleaning Sequence

Polar Samples	Non-polar Samples**	Proteinaceous Samples
1. water	1. isopropanol (or an appropriate isopropanol/ water mixture*)	Option 1: Inject repeated aliquots of dimethylsulfoxide (DMSO)
2. methanol	2. tetrahydrofuran (THF)	Option 2: gradient of 10% to 90% B where: A = 0.1% trifluoroacetic acid (TFA) in water B = 0.1% trifluoroacetic acid (TFA) in acetonitrile (CH ₃ CN)
3. tetrahydrofuran (THF)	3. dichloromethane	
4. methanol	4. hexane	
5. water	5. isopropanol (followed by an appropriate isopropanol/ water mixture*)	Option 3: Flush column with 7M guanidine hydrochloride, or 7M urea
6. mobile phase	6. mobile phase	

* Use low organic solvent content to avoid precipitating buffers.

** Unless a Hexane Tetrahydrofuran Compatibility Kit (PN 205000464) has been installed, running solvents such as THF or hexane should only be considered when the column cannot be cleaned by running neat, reversed-phase organic solvents such as acetonitrile. Reduce flow rate, lower operating temperatures and limit system exposure to THF and/or hexane.

b. For Normal-Phase Conditions

The HSS Cyano column can be used for both reversed-phase separations as well as normal-phase separations. The column is originally shipped in acetonitrile and is ready to use for reversed-phase conditions.

If you intend to use the column for normal-phase applications, you will need to condition the column with the following procedure:

1. Flush the column with a minimum of 20 column volumes of 100% methanol using a low flow rate to avoid overpressuring the LC system. Refer to Table 1 for minimum solvent volume.
2. Flush the column with a minimum of 20 column volumes of 100% isopropanol using a low flow rate to avoid overpressuring the LC system. Refer to Table 1 for the minimum solvent volume.
3. Flush the column with a minimum of 20 column volumes of 100% dichloromethane using a low flow rate to avoid overpressuring the LC system. Refer to Table 1 for the minimum solvent volume.
4. Flush the column with the intended mobile-phase conditions until a stable baseline is achieved.

c. Storage

For periods longer than four days at room temperature, store reversed-phase ACQUITY UPLC HSS columns in 100% acetonitrile. For elevated temperature applications, store immediately after use in 100% acetonitrile for the best column lifetime. Do not store columns in buffered eluents. If the mobile phase contained a buffer salt, flush the column with 10 column volumes of HPLC grade water (see Table 1 for common column volumes) and replace with 100% acetonitrile for storage. Failure to perform this intermediate step could result in precipitation of the buffer salt in the column when 100% acetonitrile is introduced. Completely seal column to avoid evaporation and drying out of the bed.

Note: If a column has been run with a mobile phase that contains formate (e.g., ammonium formate, formic acid, etc.) and is then flushed with 100% acetonitrile, slightly longer equilibration times may be necessary when the column is re-installed and run again with a formate-containing mobile phase.

d. Storage post Normal-Phase use

For rapid equilibration upon start-up, it is recommended that you store the HSS Cyano column in the mobile phase that is commonly used for your normal-phase separation. Completely seal the column to avoid evaporation and drying out of the bed.

IV. INTRODUCING ECORD INTELLIGENT CHIP TECHNOLOGY

a. Introduction

The eCord intelligent chip is a new technology that will provide the history of a column's performance throughout its lifetime. The eCord will be permanently attached to the column to assure that the column's performance history is maintained in the event that the column is moved from one instrument to another.

At the time of manufacture, tracking and quality control information will be downloaded to the eCord. Storing this information on the chip will eliminate the need for a paper Certificate of Analysis. Once the user installs the column, the software will automatically download key parameters into a column history file stored on the chip. In this manual, we explain how the eCord will provide a solution for easily tracking the history of the columns, reduce the frustration of paperwork trails, and give customers the reassurance that a well-performing column is installed onto their instruments.

b. Installation

Install the column into the column heater. Plug the eCord into the side of the column heater. Once the eCord is inserted into the column heater the identification and overall column usage information will be available allowing the user to access column information on their desktop.

c. Manufacturing Information

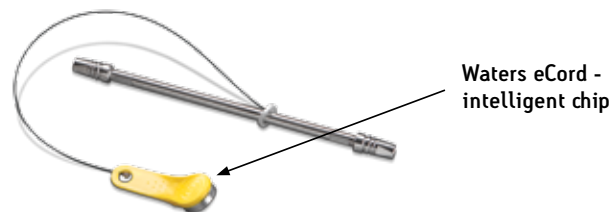


Figure 1. Waters eCord Intelligent Chip


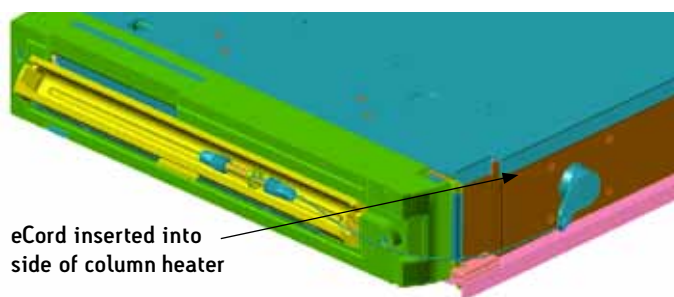


Figure 2. eCord Inserted into Side of Column Heater



The eCord chip provides the user with an overview of the bulk material QC test results.

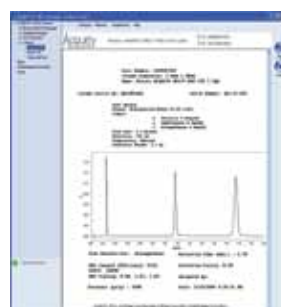

The eCord chip provides the user with QC test conditions and results on the column run by the manufacturer. The information includes mobile phases, running conditions and analytes used to test the columns. In addition the QC results and acceptance is placed onto the column.

d. Column Use Information

The eCord chip provides the customer with column use data. The top of the screen identifies the column including chemistry type, column dimensions and serial number. The overall column usage information includes the total number of samples, total number of injections, total sample sets, date of first injection, date of last injection, maximum pressure and temperature. The information also details the column history by sample set including date started, sample set name, user name, system name, number of injections in the sample set, number of samples in the sample set, maximum pressure and temperature in the sample set and if the column



V. ADDITIONAL INFORMATION

a. Tips for Maximizing ACQUITY UPLC HSS Column Lifetimes

- To maximize ACQUITY UPLC column lifetime, pay close attention to:
 - Water quality (including water purification system)
 - Solvent quality
 - Mobile phase preparation, storage and age
 - Sample, buffer and mobile phase solubilities
 - Sample quality and preparation.
- When problems arise, often only one improper practice must be changed.
- Always remember to:
 - Use in-line filter unit or, preferably, a VanGuard Pre-Column.
 - Discourage bacterial growth by minimizing the use of 100% aqueous mobile phases where possible.
 - Change aqueous mobile phase every 24 – 48 hours (if 100% aqueous mobile phase use is required).
 - Discard old 100% aqueous mobile phases every 24-48 hours to discourage bacterial growth.
 - Add 5% - 10% organic modifier to mobile phase A and adjust gradient profile.
 - Filter aqueous portions of mobile phase through 0.2 µm filter.
 - Maintain your water purification system so that it is in good working order.
 - Only use ultra pure water (18 megohm-cm) water and highest quality solvents possible. HPLC grade water is not UPLC grade water.
- Consider sample preparation (e.g., solid-phase extraction, filtration, etc).
- Avoid (where possible):
 - 100% aqueous mobile phases (if possible)
 - HPLC-grade bottled water
 - “Topping off” your mobile phases
 - Old aqueous mobile phases. Remember to rinse bottles thoroughly and prepare fresh every 24 to 48 hrs
 - Using phosphate salt buffer in combination with high ACN concentrations (e.g., > 70%) due to precipitation.
- Don't: assume a “bad” column is the culprit when high backpressure or split peaks are observed:
 - Investigate cause of column failure
 - Backpressure
 - Mobile phase(s), bacteria, precipitation and/or samples
 - Peak splitting
 - Sample quality
 - Injection solvent strength.
- Remember: the diameter of UPLC columns (1.0, 2.1 and 3.0 mm ID) are often lower than that of a conventional HPLC column and therefore, mobile phases last much longer. To reduce the chances of mobile-phase contamination or degradation, only prepare what you need for analysis or store excess bulk quantities in a refrigerated environment.

7. Mobile phase-related questions to ask:

- Am I using 100% aqueous mobile phases? Am I able to add a small amount of organic modifier to my mobile phase A?
- Do I filter my aqueous mobile phases through 0.2 µm filters?
- How old is my mobile phase? Do I label the bottle with preparation date?
- Do I “top off” or do I prepare fresh mobile phases every 24 – 48 hrs?
- What is the quality of my water? Has the quality recently changed? How is my water purification system working? When was it last serviced?

- Am I working with pH 7 phosphate buffer (which is VERY susceptible to bacterial growth)?

8. Sample-related questions to ask:

- If I inject neat standards prepared in mobile phase do I observe these problems?
- If I prepare my standards in water and prepare them like samples (e.g. SPE, filtration, etc.) do I still observe these problems?
- Has the quality of my samples changed over time?

b. Recommended Flow Rates and Backpressures for Reversed-Phase ACQUITY UPLC HSS Columns

1.0 mm i.d. Columns (40 °C)								
UPLC Linear Velocity (mm/sec)	3		4		5		6	
Column Dimensions	Flow Rate (mL/min)	Backpressure (psi)	Flow Rate (mL/min)	Backpressure (psi)	Flow Rate (mL/min)	Backpressure (psi)	Flow Rate (mL/min)	Backpressure (psi)
1.0 x 50 mm	0.1	3800	0.13	5600	0.17	6500	0.2	7600
1.0 x 100 mm	0.1	7500	0.13	9800	0.17	12800	0.2	15100
1.0 x 150 mm	0.1	11300	0.13	14600	0.17	19100	0.2	25500

2.1 mm i.d. Columns (40 °C)								
UPLC Linear Velocity (mm/sec)	3		4		5		6	
Column Dimensions	Flow Rate (mL/min)	Backpressure (psi)	Flow Rate (mL/min)	Backpressure (psi)	Flow Rate (mL/min)	Backpressure (psi)	Flow Rate (mL/min)	Backpressure (psi)
2.1 x 30 mm	0.45	2700	0.60	3600	0.75	4500	0.9	5400
2.1 x 50 mm	0.45	4200	0.60	5700	0.75	7100	0.9	5800
2.1 x 100 mm	0.45	8000	0.60	10700	0.75	13400	0.9	16100
2.1 x 150 mm	0.45	11800	0.60	15800	0.75	19700	0.9	23700

Note: 1) ACQUITY HSS UPLC 1.8 µm particle reversed-phase columns
 2) ACN/Aqueous gradient, Pmax at ~30% ACN
 3) Approximate maximum total system backpressure given.

3.0 mm i.d. Columns (40 °C)								
UPLC Linear Velocity (mm/sec)	3		4		5		6	
Column Dimensions	Flow Rate (mL/min)	Backpressure (psi)	Flow Rate (mL/min)	Backpressure (psi)	Flow Rate (mL/min)	Backpressure (psi)	Flow Rate (mL/min)	Backpressure (psi)
3.0 x 30 mm	0.9	3100	1.17	4000	1.53	5300	1.8	6200
3.0 x 50 mm	0.9	4600	1.17	6000	1.53	7800	1.8	9200
3.0 x 100 mm	0.9	8300	1.17	10800	1.53	14100	1.8	16700
3.0 x 150 mm	0.9	12000	1.17	15700	1.53	20500	1.8	24100

Sales Offices

**Austria and European Export
(Central South Eastern Europe, CIS
and Middle East)** 43 1 877 18 07

Australia 61 2 9933 1777

Belgium 32 2 726 1000

Brazil 55 11 4134 3788

Canada 1 800 252 4752

China 86 10 5293 6688

Czech Republic 420 2 617 11384

Denmark 45 46 59 8080

Finland 358 9 5659 6288

France 33 1 30 48 72 00

Germany 49 6196 400 600

Hong Kong 852 2964 1800

Hungary 36 1 350 5086

India 91 80 2837 1900

Ireland 353 1 448 1500

Italy 39 02 265 0983

Japan 81 3 3471 7191

Korea 82 2 6300 4800

Mexico 52 55 52 00 1860

The Netherlands 31 76 508 7200

Norway 47 6 384 6050

Poland 48 22 833 4400

Puerto Rico 1 787 747 8445

Russia/CIS 7 495 727 4490/ 290 9737

Singapore 65 6593 7100

Spain 34 93 600 9300

Sweden 46 8 555 115 00

Switzerland 41 56 676 7000

Taiwan 886 2 2501 9928

United Kingdom 44 208 238 6100

All other countries:

Waters Corporation U.S.A.

1 508 478 2000

1 800 252 4752

www.waters.com

Waters

THE SCIENCE OF WHAT'S POSSIBLE.™

©2009 Waters Corporation. Waters, The Science of What's Possible, ACQUITY, ACQUITY UPLC, UPLC, eCord, VanGuard, Oasis and Sep-Pak are trademarks of Waters Corporation. All other trademarks are the property of their respective owners.

October 2011 715001429 Rev C KK-PDF

Waters Corporation
34 Maple Street
Milford, MA 01757 U.S.A.
T: 1 508 478 2000
F: 1 508 872 1990
www.waters.com

